

FORMULATED CHEMISTRIES FOR PROCESS CLEANING

PHARMACEUTICAL DETERGENTS

Developing a validated cleaning process for a large stainless steel bioreactor or for tablet manufacturing equipment will present unique challenges. STERIS offers a complete line of pharmaceutical detergents for the most challenging validated clean-in-place (CIP), clean-out-of-place (COP), and manual cleaning applications. You can be assured that our detergents will deliver consistent cleaning performance so there will be no surprises in your validated manufacturing process.

Common Applications

STERIS pharmaceutical detergents deliver unmatched performance in the following critical applications:

- CIP cleaning of large vessels and process piping.
- COP cleaning in pharmaceutical washers, recirculation baths and ultrasonic tanks.
- Manual Cleaning of complex equipment with sensitive substrates.
- Rouge removal and passivation of stainless steel process equipment.

Technical Support

Our technical support services include the Process And Cleaner Evaluation (PACE®) program, which is an evaluation service designed to provide our Customers with recommendations for an effective cleaning process. Once an evaluation has been completed, STERIS provides a report that assists Customers in developing a cleaning process by defining parameters based on chemistry, concentration, cleaning time, temperature, cleaning method and water quality. The PACE program is an essential first step for any cleaning application.

A highly qualified, industry-recognized team of chemists, microbiologists and engineers is available to offer product and process consultation. STERIS Technical Support currently provides both on- and off-site seminars with topics focusing on process cleaning and cleaning validation. An extensive library of technical data, laboratory reports, analytical methods and case studies have been developed including toxicity and substrate compatibility studies.



Products

STERIS pharmaceutical detergents are designed for your most difficult cleaning challenges. They employ multiple cleaning mechanisms such as wetting, solubility, dispersion, emulsification, hydrolysis and oxidation which work together to penetrate, dissolve and rinse away process residues. Each product is supported by our industry leading technical support staff as well as an extensive documentation package to help facilitate cleaning validation.

Alkaline Detergents

Alkaline Detergents are used for the majority of cleaning applications due to their ability to penetrate and efficiently remove a wide range of organic residues.

- CIP 100® Alkaline Detergent
 - Formulated for a broad range of residues and applications
- CIP 130™ Alkaline Detergent
 - Formulated with a multiple surfactant system
- CIP 150® Chlorinated Alkaline Detergent
 - Formulated with an oxidizing agent
- CIP Neutralizer™ Alkaline Based Neutralizer
 - Increases the pH of a solution
- DA-7645® Alkaline Detergent
 - Formulated for manual cleaning
- Foam 140® Alkaline Foaming Detergent
 - Produces heavy foam with spray applicator
- ProKlenz® ONE Alkaline Detergent
 - Formulated for a broad range of residues and applications
 - EPA registered disinfectant and virucide

Neutral Detergents

Neutral Detergents are commonly used in manual cleaning applications or with sensitive substrates due to their mild pH and enhanced compatibility. They may also be used to enhance the performance of other detergents.


- CIP 300® Neutral Detergent
 - Surfactant based formulation
- CIP Additive™ Detergent
 - Surfactant based formulation
- ProKlenz Booster Detergent
 - Surfactant based formulation with oxidizing agent
- ProKlenz NpH Neutral Detergent
 - Surfactant and solvent based formulation



Acid Detergents

Acid Detergents are commonly used for the removal of inorganic residues and salts. Additionally, they are widely used for rouge removal and passivation of stainless steel equipment.

- CIP 200® Acid Detergent
 - Phosphoric acid formulation
- CIP 220® Acid Detergent
 - Organic acid formulation
- Foam 240® Acid Foaming Detergent
 - Produces heavy foam with spray applicator
- ProKlenz TWO Acid Detergent
 - Organic acid formulation
 - Formulated for rouge removal and passivation

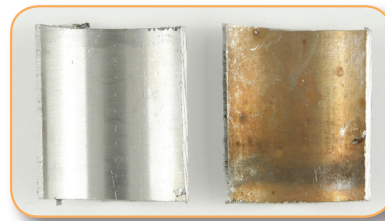

STERIS has a comprehensive offering of detergents, disinfectants, skin care products and sterility assurance products that support your needs. We also have world class technical support to design the most effective cleaning program for your facility. Contact your local sales representative or refer to our website at www.sterislifesciences.com

 **STERIS**®
STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 USA
www.sterislifesciences.com

STERIS OFFICES WORLDWIDE

Asia Pacific	+65 (0) 68 41 7677	Latin America	+1 800 444 9009
Canada	+1 800 661 3937	Spain	+34 (0) 916 585 920
France	+33 (0) 2 38 70 83 50	United Kingdom	+44 (0) 1256 840400
Germany	+49 (0) 221 466120	United States	+1 800 444 9009
Italy	+39 (0) 2 21303 424		

Document # M3488EN
Printed 01/2011

©2011 STERIS Corporation.
All rights reserved.